



*available options shown

- FAST CHAIR PREMIUM can be used as a stair chair, evacuation chair or transport chair
- Stair-track easily attachable or removable to minimize carried by the operator weight
- Can be also used with powered stair-track Powertraxx (optional accessory)
- Lightweight all-aluminum frame supports patients up to 250 kg
- Fold-out comfortable armrests and side handle protective guards
- Smooth and quiet wheels - two 125 mm front swivel and two 150 mm rear fixed wheels with locks
- Head support, three restraints, armrests and comfort footbar included
- Wide rear locking handles that are angled for better ergonomomy
- Plastic seat and backrest are easy to clean

For all product variations please see the price list

Item Number	60-4150-001
Height (min)	1000 mm / 39.37 in
Width (overall)	585 mm / 23.03 in
Width (seat)	420 mm / 16.54 in
Length (folded)	1000 mm / 39.37 in
Width (folded)	585 mm / 23.03 in
Depth (folded)	240 mm / 9.45 in
Weight (w/o track)	13.8 kg / 30.42 lbs
Weight (with track)	20.7 kg / 45.64 lbs
Load Limit	250 kg / 551.16 lbs
Power Supply	28V LI-Ion Milwaukee Battery



Lockable compact folding for easy transport or storage with easy locking and unlocking



Stair track can be easily attached to or removed from the chair (to minimize the the weight that needs to be carried by the operator)



Wheels have minimum debris collection while maintaining good traction and comfort, and provide very good maneuverability in confined spaces such as stair landings or tight corridors



Top telescoping handle is extendable and adjustable in four positions (1000 to 1580 mm)



Front telescopic handles adjustable in five positions (maximum extension of 340 mm)

OPTIONS

FAST/FST WALL MOUNT - designed to securely hold the FAST Chair in both a static or a mobile application. Single hand release allows for quick chair fastening or detachment.

60-4190-001




4-POINT RESTRAINT for enhanced safety

62-4400-001



BIG WHEELS FOR FAST CHAIR can be quickly and easily attached or removed and stored

60-0407-001

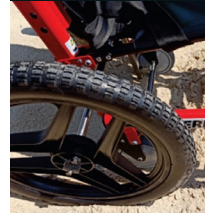


FAST POWERTRAXX intelligent battery powered stair-track system that glides up to 250kg loads up or down the stairs without carrying or lifting.

60-0403-001




MESH BASKET for storing accessories such as oxygen bottle

65-4123-001

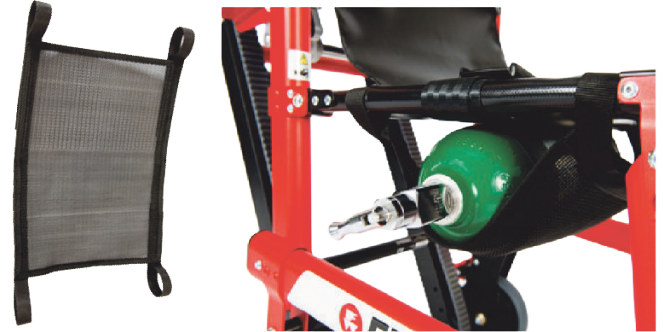

FAST-LIFT

- suitable for FAST PREMIUM, FAST, and FST chairs
- provides spring-assisted lifting (for example onto ambulance wall or doors)
- safely secures the chair

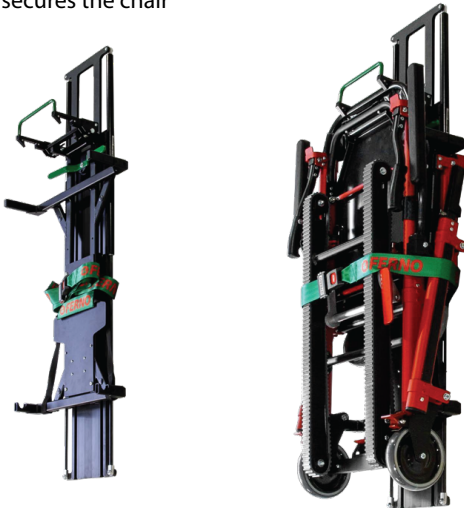

60-0254-002

FAST-LIFT X

- identical to FAST-LIFT
- with adjustable forward/backward extension for recessed storage spaces
- safely secures the chair



60-0254-001

FAST CHAIR PREMIUM complies with and exceeds the requirements of the international standard EN 1865.

